


# Vision Tester VT-10



TOPCON CORPORATION

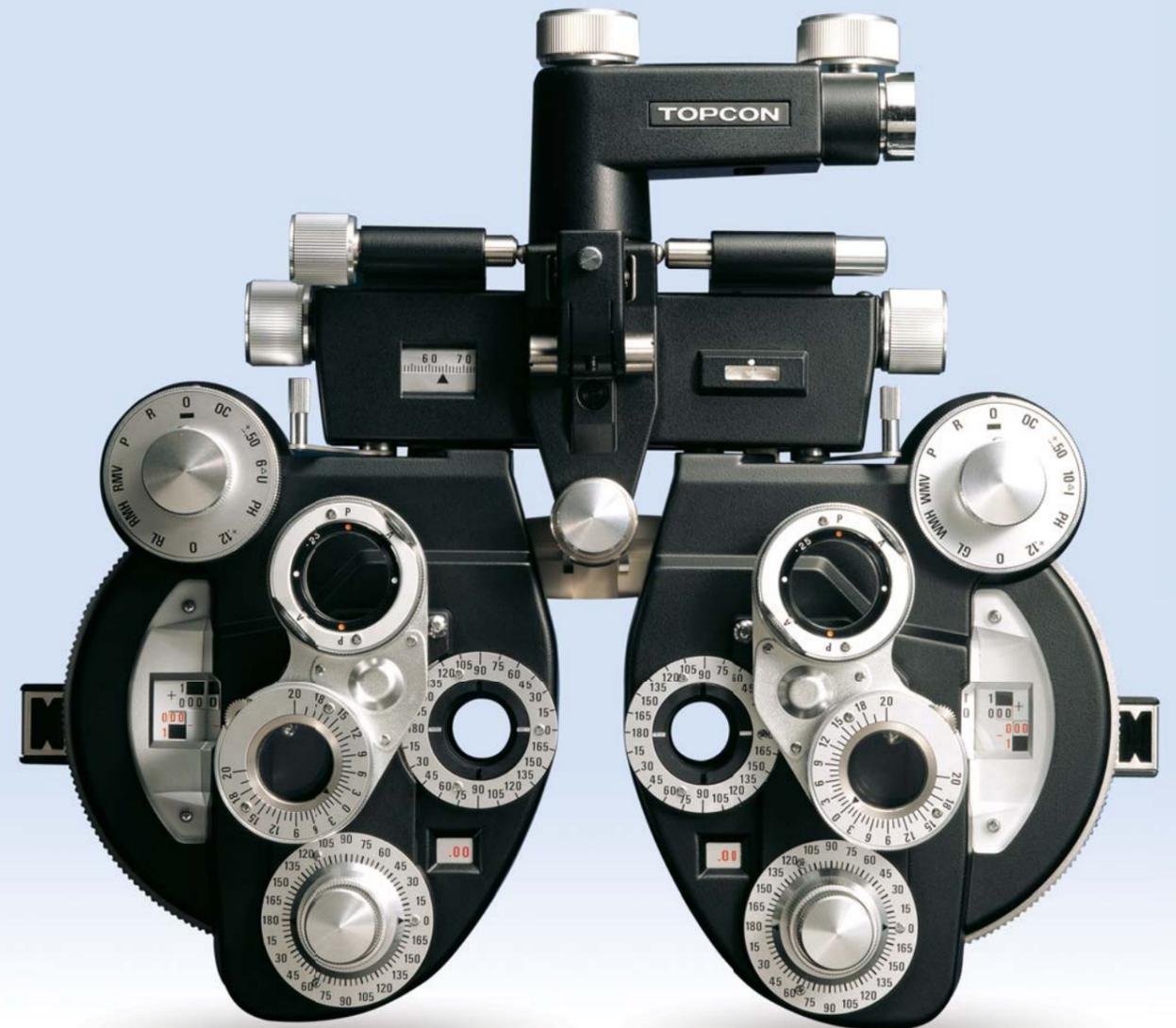

TOPCON CORPORATION



EC9711081

\*Subject to change in design and/or specifications without advanced notice.  
In compliance with the terms of the Export Administration Regulation of the United States of America, this product may not be available in some regions or countries.

**IMPORTANT** In order to obtain the best results with this instrument, please be sure to review all user instructions prior to operation.

**TOPCON CORPORATION**  
CONNECTING VISIONS 75-1 Hasunuma-cho, Itabashi-ku, Tokyo 174-8590, Japan. Phone: 3-3558-2523/2522 Fax: 3-3965-6898 www.topcon.co.jp

**TOPCON MEDICAL SYSTEMS, INC.**  
111 Bauer Drive, Oakland, NJ 07436, U.S.A.  
Phone: +1-201-599-5100 Fax: +1-201-599-5250 www.topconmedical.com

**TOPCON CANADA INC.**  
110 Provencher Avenue, Boisbriand, QC J7G 1N1 CANADA  
Phone: +1-450-430-7771 Fax: +1-450-430-6457 www.topcon.ca

**TOPCON EUROPE MEDICAL B.V.**  
Esebaan 11, 2908 LJ Capelle a/d IJssel, P.O.Box 145,  
2900 AC Capelle a/d IJssel, THE NETHERLANDS  
Phone: +31-(0)10-4585077 Fax: +31-(0)10-4585045  
E-mail: medical@topcon.nl; www.topcon.eu

**ITALY OFFICE**  
Viale dell'Industria 60, 20037 Paderno Dugnano, (Milano), ITALY  
Phone: +39-039-9186871 Fax: +39-02-91081091 E-mail: topconitaly@tiscali.it; www.topcon.it

**DENMARK OFFICE**  
Præstemarksvej 25, 4000 Roskilde, DENMARK  
Phone: +45-46-327500 Fax: +45-46-327555  
E-mail: topcon@topcondenmark.dk www.topcondenmark.dk

**IRELAND OFFICE**  
Unit 276, Blanchardstown, Corporate Park 2 Ballycoolin Dublin 15, IRELAND  
Phone: +353-18975900 Fax: +353-18293915 E-mail: medical@topcon.ie; www.topcon.ie

**TOPCON S.A.R.L.**  
89, rue de Paris 92585 Clichy, FRANCE  
Phone: +33(0)1-41069494 Fax: +33(0)1-47390251 E-mail: topcon@topcon.fr; www.topcon.fr

**TOPCON DEUTSCHLAND G.m.b.H.**  
Hanns-Martin-Schleyer Strasse 41, D-47877 Willich, GERMANY  
Phone: +49-(0)2154-8850 Fax: +49-(0)2154-885177 E-mail: medical@topcon.de; www.topcon.de

**TOPCON SCANDINAVIA A.B.**  
Neogatan 2, P.O.Box 25, 43151 Mölndal, SWEDEN  
Phone: +46-(0)31-7109200 Fax: +46-(0)31-7109249 E-mail: medical@topcon.se; www.topcon.se

**TOPCON ESPAÑA S.A.**  
**HEAD OFFICE**  
Frederic Mompou 4 Esc. A Bajos 3, 08960 Sant Just Desvern Barcelona, SPAIN  
Phone: +34-935-4734057 Fax: +34-93-4733332 E-mail: medical@topcon.es; www.topcon.es  
**PORTUGAL OFFICE**  
Rua da Forte, 6-6A, L-0, 22, 2790-072 Carnaxide, PORTUGAL  
Phone: +351-210-994626 Fax: +351-210-938786 www.topcon.pt

**TOPCON ( GREAT BRITAIN ) LTD.**  
Topcon House, Kennet Side, Bone Lane, Newbury, Berkshire RG14 5PX United Kingdom  
Phone: +44-(0)1353-551120 Fax: +44-(0)1353-551170 E-mail: medical@topcon.co.uk; www.topcon.co.uk

**TOPCON POLSKA Sp. z o.o.**  
ul. Warszawska 23, 02-470 Siewierz, POLAND  
Phone: +48-(0)32-6705045 Fax: +48-(0)32-6713405 www.topcon-polska.pl

**TOPCON SINGAPORE MEDICAL PTE. LTD.**  
60 Alexandra Terrace, #08-27 The Comtech, Singapore 118502  
Phone: +65-68726060 Fax: +65-67736150 E-mail: medical\_sales@topcon.com.sg www.topcon.com.sg

**TOPCON INSTRUMENTS ( MALAYSIA ) SDN.BHD.**  
No. D1, (Ground Floor), Jalan Excela 2, Off Jalan Ampang Putra,  
Taman Ampang Hilir, 55100 Kuala Lumpur, MALAYSIA  
Phone: +60-(0)3-42709866 Fax: +60-(0)3-42709766

**TOPCON INSTRUMENTS ( THAILAND ) CO.,LTD.**  
77/162 Simsathorn Tower, 37th Floor, Krungthorbnuri Rd., Klongtongnai,  
Klongsarn, Bangkok 10600, THAILAND  
Phone: +66(0)2-440-1152-7 Fax: +66(0)2-440-1158

**TOPCON CORPORATION BEIJING OFFICE**  
Block No.9, Kangding Street Beijing Economic-Technological Development Area,  
Beijing, 100176, CHINA  
Phone: +86-(0)10-6780-2799 Fax: +86-(0)10-6780-2790

**TOPCON CORPORATION SHANGHAI OFFICE**  
14L Huamin Empire Plaza, No.726, Yan'an Xi Road,  
Shanghai, 200050, CHINA  
Phone: +86-(0)21-5238-7722 Fax: +86-(0)21-5237-0761

**TOPCON CORPORATION BEIRUT OFFICE**  
P.O.Box 70-1002 Antelias, Beirut, LEBANON  
Phone: +961-4-523525/523526 Fax: +961-4-521119

**TOPCON CORPORATION DUBAI OFFICE**  
P.O.Box 293705, Dubai Airport Free Zone L.L.U., J-12, Dubai, U.A.E  
Phone: +971-4-299-5900 Fax: +971-4-299-5901

## Versatility & Convenience Combined with Unparalleled Response & Precision

Precision and convenience are a tradition in high quality TOPCON optical instruments. The VT-10 vision tester advances this tradition still further, making refraction quicker and easier than ever before, providing unequaled performance and exceptional versatility. Precision in design and workmanship guarantees accuracy in refraction measurements—accuracy that won't deteriorate even after long years of use. Total quality from TOPCON means you can have total confidence in the performance you'll enjoy from the VT-10 vision tester.



**VT-10**  
Vision Tester

## A Half-Century of Dedication to Superior Optical Instrumentation

TOPCON has remained a leader and major innovator in optical instrumentation and technology for more than half a century. Intensive dedication to engineering research and manufacturing precision has made TOPCON the name to know for reliable performance.

### Precision Engineering for Superior Performance;

Designed for total precision, the VT-10 vision tester represents TOPCON optical engineering at its very best—extreme accuracy over the widest possible testing range. Internal mechanisms are precisely machined to the finest tolerances; parts and components are made from the finest materials, all corrosion-resistant and tested for performance and durability.

### Multi-Coated Lenses for Superior Optical Performance

Every lens used in the VT-10 receives a special multi-coating that eliminates flare, ghost images and reflections. In addition, the multi-coating assures clearer, brighter images with improved contrast and excellent color fidelity.



The TOPCON VT-10 is available in either black or gray finishes.

# Controls & Scales Designed to Make Refraction Quick & Comfortable; Features for Functional Versatility



## Wide Testing Range

The testing range available with the VT-10 is extraordinary, the standard cylinder power range extends to 6.00D in 0.25D increments. In addition, auxiliary lens sets are incorporated to fine-tune cylinder power measurements to 0.12D increments or extend cylinder power range to 8.00D. Accurate prism reading even below  $0.5\Delta$  are possible thanks to the  $20\Delta$  range ( $1\Delta$  increments) of the precision-gear rotary prism loupes. The VT-10 is also available in plus or minus cylinder models.




## Synchronized Cross-Cylinders

The most significant time-saving feature of the VT-10 is its synchronization of cross-cylinder loupe rotation to cylinder axis adjustment. Special gearing automatically synchronizes the rotation of the loupe each time you change cylinder axis correction. The result? You simply read a single scale, one time, after completing a full series of axis checks. The bi-loupe turret lets you index power or axis loupes with a single touch.



## Easy Controls & Clear Scales

Control of cylinder axis and power testing is amazingly simple; concentric control knobs within a double-calibration axis scale make adjustment a quick, one hand action through the full  $360^\circ$  range. The axis scale reads in  $5^\circ$  increments, with  $15^\circ$  markings. An additional axis scale is located at the main aperture, for convenient reference during retinoscopy. The low position of the control knobs keeps the operator's hands out of the line of vision.

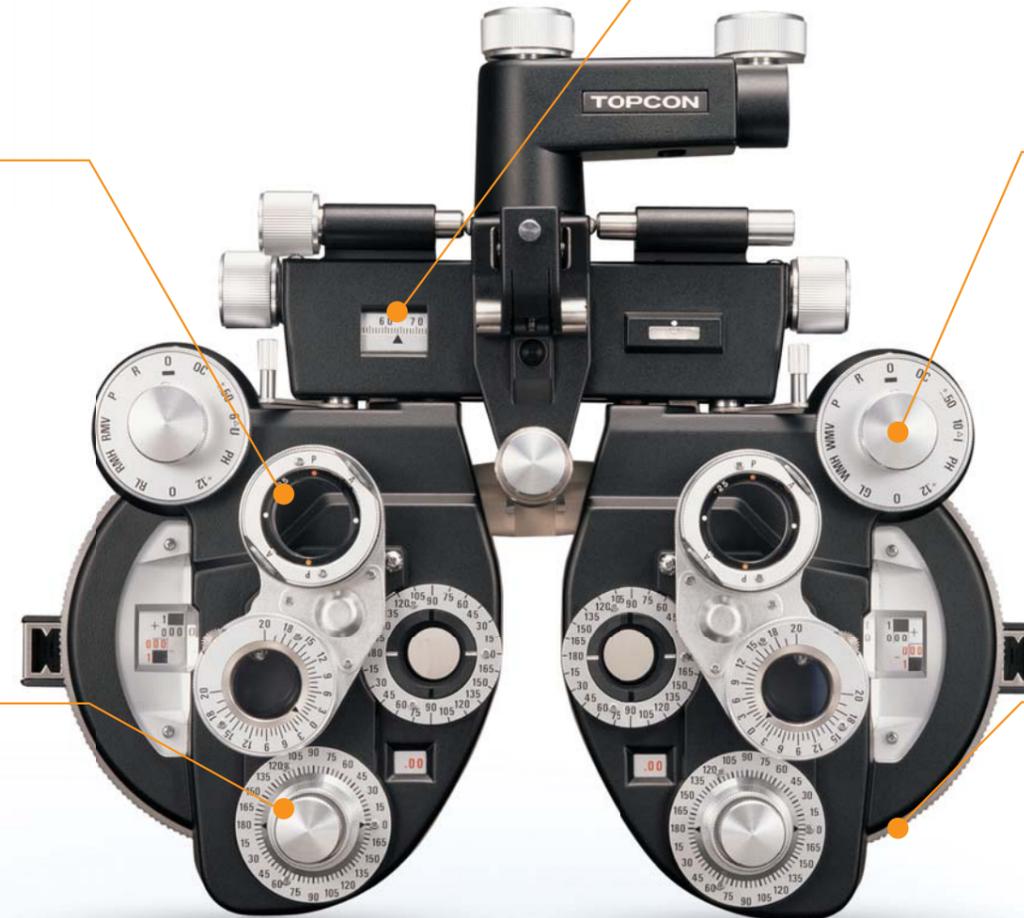

## Versatile Adjustments

Convenient dual controls make it easy for the operator to adjust interpupillary distance, in a range from 48mm to 75mm. The mechanism is sealed for protection and smooth functioning. An easy-to-read scale provides 1mm graduations, and a vergence adjustment controls reductions in P.D. during near testing.



## A Wide Range of Versatile Auxiliary Lenses

The TOPCON VT-10 comes with a full complement of versatile auxiliary lenses—10 different lenses plus two open apertures, any of which can be indexed into position with one-touch speed and convenience. A red lens/green lens shift, for example, is immediately available for binocular vision testing. The auxiliary lens plate has a flat black finish to prevent stray reflections during use.



## Fast, Precise Results

TOPCON VT-10 incorporates a single-dial power control that lets the operator quickly dial in any spherical power from +16.75D through 0 to -19.00D, in 0.25D increments. To make settings faster, a strong sphere control allows adjustment in 3.00D steps; the auxiliary lens knob provides for fine-tuning in 0.12D increments for absolute accuracy.




# Other Features

## Simple Convergence System

Precise measurements at both near and distant testing for optimal results.



One quick adjustment of the dual convergence levers and the VT-10 provide total accuracy in near vision testing measurement.

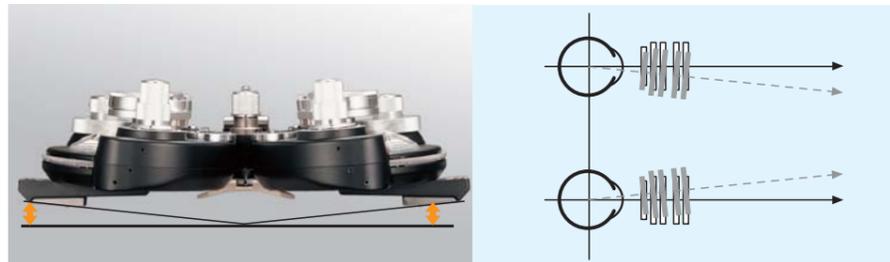

## Sealed Optics

The examination apertures feature cover lenses for convenience and protection. Dust, eyelashes, etc. can't get in, and the lenses are easily removed for simple cleaning. The patient always has a bright, clear image to view. (Standard accessories include face shields that can be sterilized for optimum comfort and hygiene.)



## Accessories

The holder for rotary near point card mounts to a sturdy rod with a distance scale from 5 to 28 inches—centimeter and diopter scales are included. A hinged bracket moves it into position immediately, or raises it out of the way. Accessory lenses come in an accessory box, holding dual sets or -0.12D and -2.00D cylinder lenses.



## Specifications

**Type** Subjective refraction testing equipment

**Item and range of measurement** The value in ( ) is with accessory lens or auxiliary lens attached:

<b>Myopia</b>	0 to -19.00D in 0.25D step (0.12D step)
<b>Hyperopia</b>	0 to +16.75D in 0.25D step (0.12D step)
<b>Astigmatism</b>	0 to -6.00D ( -8.00D) in 0.25D step (0.12D step)
<b>Cross-Cylinder</b>	± 0.25D ( ± 0.37D and ± 0.50D as option)
<b>Phoria and vergency</b>	20△ to 0 to 20△ by 1△ step Inner 40△ Outer 40△ (50△) Upward and Downward 40△ (46△)
<b>Near vision test</b>	Main optical axes of both lenses converges by means of a tilt mechanism at 40cm in front of the eyes at 64mm P.D. If P.D. for far vision is less than 54mm, converging operation stops.
<b>Binocular balance test</b>	Rotary prism and polarizing lens methods are available.
<b>P.D.</b>	48 to 75mm by 1mm step
<b>Level adjustment</b>	More than 4.5° for right and left banks
<b>Dimensions</b>	335 (at 64mm P.D. setting) x 97 x 294mm
<b>Weight</b>	4.5kg

### Content of Auxiliary Lens Dial:

"O/O"	Open Aperture
"R"	Retinoscopic Lens: +1.50D standard +2.00D option
"P"	Polarizing Filter
"WMV/RMV"	Maddox Rod, vertical: white (left) red (right)
"WMH/RMH"	Maddox Rod, horizontal: white (left) red (right)
"GL"	Green Lens (left)
"RL"	Red Lens (right)
" +.12"	+0.12D Spherical
"PH"	Pin Hole
"10 I"	10 Diopter Prism Base In (left)
"6 U"	6 Diopter Prism Base Up (right)
" ± .50"	± 0.50D Fixed Cross Cylinder
"OC"	Occluder

Plus and minus cylinder models available, also the choice of color—black or light grey.

### Optional Accessories



± 0.37 / ± 0.50D cross cylinder, + 2.00D retinoscopic lens.